

Installation and User Guide

AudioVision

Capstone project

15th May, 2024



AudioVision User Manual

Table of Contents

1. Introduction
2. System Requirements
3. Installation Instructions
4. Getting Started with AudioVision
5. Frequently Asked Questions (FAQ)
6. Contacting Support



Introduction

Welcome to the user manual for AudioVision!

AudioVision is a revolutionary application designed to assist visually impaired individuals in navigating their daily activities with greater ease and safety. Using advanced technology, AudioVision offers a range of features that enable users to confidently navigate their surroundings.

What Does AudioVision Offer?

AudioVision comes equipped with six main features designed to provide an intuitive and helpful experience in various situations:

1. **Real-time Object Detection:** This feature allows the app to detect and identify objects around you in real-time, providing audio information about your surroundings.
2. **Audio Command:** Control the app using just your voice. Give voice commands, and AudioVision will respond according to your instructions.
3. **Live Tracking:** This feature lets others track your location in real-time, ensuring you are always under the watchful eye of someone you trust.
4. **Share Location:** Easily share your location with friends and family, ensuring they know where you are and can assist if needed.
5. **Route Guidance:** Get clear and detailed route guidance to reach your destination. AudioVision will inform you of every turn and direction to take.
6. **Preference Language:** Choose your preferred language for using the app, making it easier and more comfortable to navigate.



We hope AudioVision will be your trusted companion in your daily activities. This guide will help you understand how to use each feature to its fullest potential, so you can get the most out of this app.

Enjoy using AudioVision, and may this app help you lead your days with more ease and independence!



System Requirements

Before installing AudioVision, ensure your device meets the following minimum system requirements:

- **Operating System:** Android 8.0 (Oreo)
- **Hardware:**
 - Minimum 4GB of RAM
 - Minimum Quad-core processor
 - Minimum 500MB of available storage space
 - Camera (for real-time object detection)
 - Microphone (for audio commands)
- **Internet Connection:** Required for route guidance, live tracking and location sharing features.

Installation Instructions

Welcome to the installation guide for AudioVision!

Follow these simple steps to install the application on your device and start using its features to assist in your daily activities.


let's do those steps:

1. Download the App:
 - Click this link: [Click Here](#) or just scan this qr code




- Click option to download file
- Go to folder where you downloaded


- Tap App already downloaded to install the app.

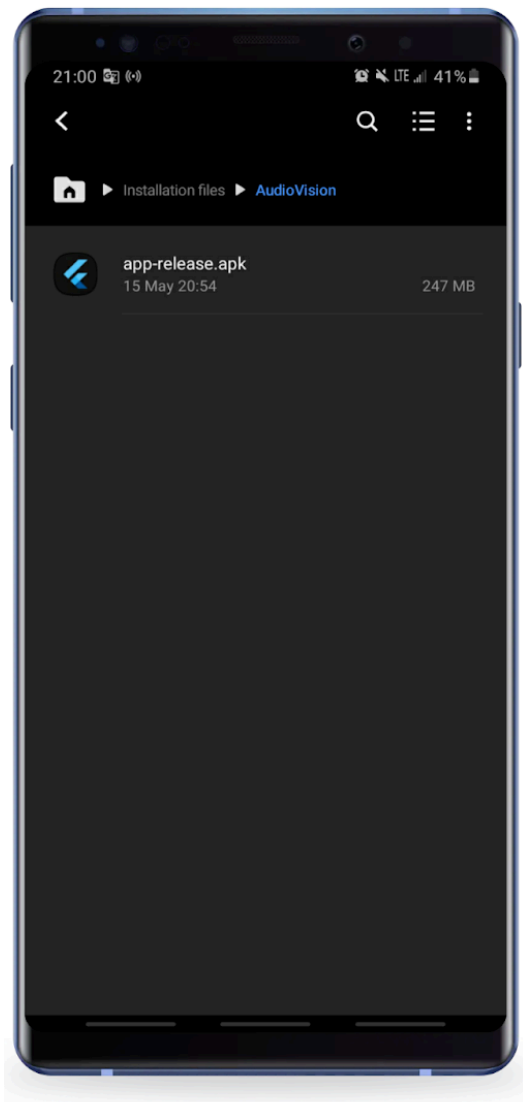

- Click "Install" to install the application

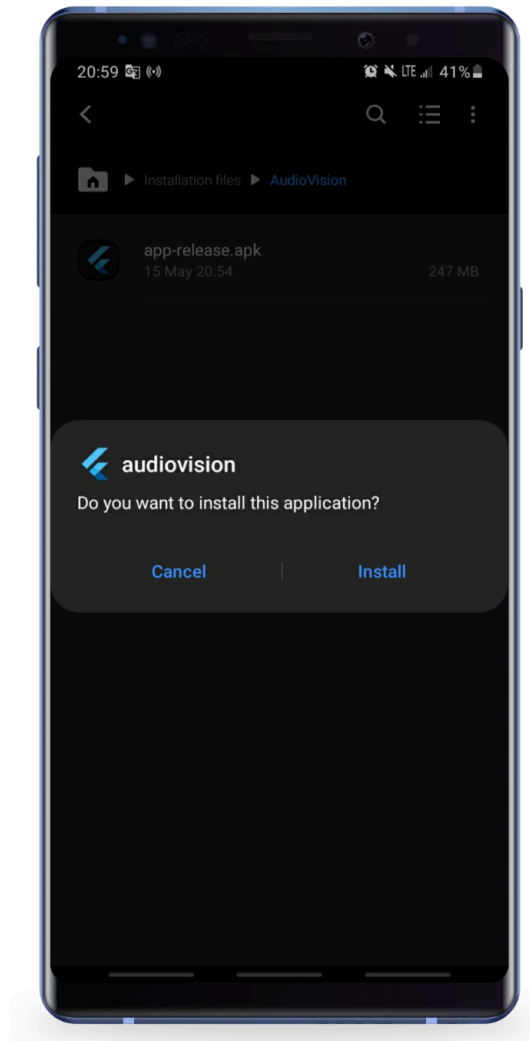

- Click "Open" to launch the AudioVision application.

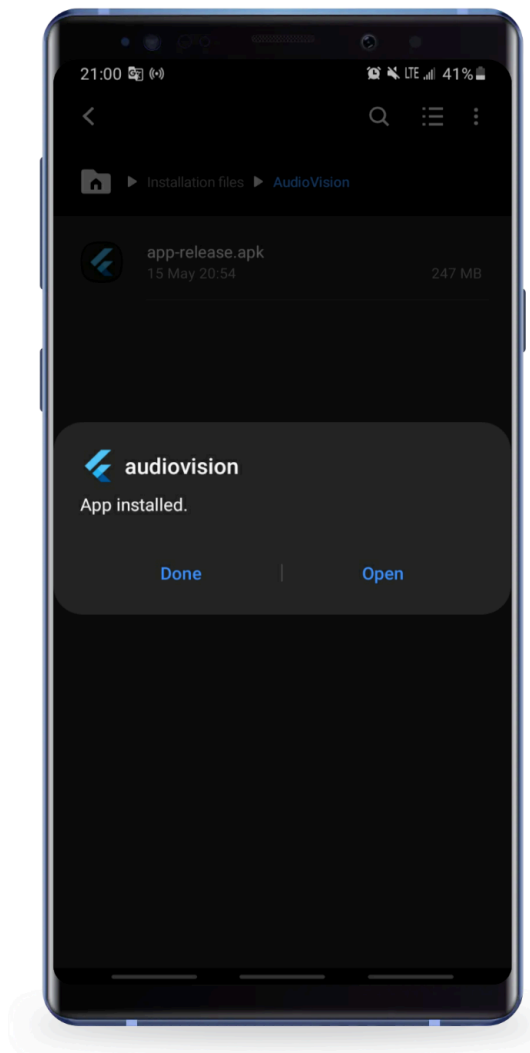



2. Grant Permissions:

- After installation, open the app.
- The app will request permissions to access your camera, microphone, and location.
- Tap "Allow" to grant these permissions.

Troubleshooting Installation Issues:

If you encounter any issues during installation, try the following solutions:

- **Insufficient Storage:** Ensure your device has enough storage space available by deleting unnecessary files or apps.
- **Internet Connection:** Verify that your device is connected to a stable internet connection.
- **Permission Denied:** Check the app permissions in your device settings and ensure all required permissions are granted.
- **Compatibility:** Ensure your device meets the minimum system requirements.

For further assistance, please contact our support team via [Our Email](#) for additional resources and help.



Getting Started with AudioVision:

Welcome to AudioVision!

Follow these steps to get started with using the app and exploring its features to assist you in your daily activities.

Voice Command

Follow these steps to get started with using the app and exploring its features to assist

Welcome to the Navigation Application User Guide! This guide will help you understand and use the various voice commands available in both English and Indonesian to interact with the app. Follow the instructions below to make the most out of your navigation experience.

English Commands

1. Setting a Destination*

- Command : go (destination) or going (destination)
- Usage : To set your desired destination.
- Example : "Go to Central Park" or "Going to Central Park"

2. Stopping Navigation

- Command : stop or exit
- Usage : To stop the current navigation.
- Example : "Stop navigation" or "Exit navigation"



3. Starting Navigation

- Command : start or navigate
- Usage : To start navigation to the set destination.
- Example : "Start navigation" or "Navigate to my destination"

4. Checking Obstacles in Front

- Command : in front
- Usage : To check for obstacles in front of you.
- Example : "What is in front of me?"

5. Getting the Next Step

- Command : next or step
- Usage : To get the next step in your navigation.
- Example : "What's the next step?" or "Next step"

6. Checking the Current Destination

- Command : destination or location
- Usage : To check your current destination.
- Example : "What is my destination?" or "Current location"

7. Sharing Location

- Command : share
- Usage : To share your location with others.
- Example : "Share my location"



8. Switching Languages

- Command : language or switch
- Usage : To switch between English and Indonesian.
- Example : "Switch language" or "Change to Indonesian"

Indonesian Commands

1. Menentukan Tujuan

- Perintah : pergi (destinasi)
- Penggunaan : Untuk menentukan tujuan yang diinginkan.
- Contoh : "Pergi ke Monas"

2. Menghentikan Navigasi

- Perintah : berhenti atau keluar
- Penggunaan : Untuk menghentikan navigasi saat ini.
- Contoh : "Berhenti navigasi" atau "Keluar navigasi"

3. Memulai Navigasi

- Perintah : mulai atau navigasi
- Penggunaan : Untuk memulai navigasi ke tujuan yang telah ditetapkan.
- Contoh : "Mulai navigasi" atau "Navigasi ke tujuan saya"

4. Memeriksa Rintangan di Depan

- Perintah : depan



- Penggunaan : Untuk memeriksa rintangan di depan Anda.

- Contoh : "Apa yang ada di depan saya?"

5. Mendapatkan Langkah Berikutnya

- Perintah : selanjutnya atau langkah

- Penggunaan : Untuk mendapatkan langkah berikutnya dalam navigasi.

- Contoh : "Langkah berikutnya apa?" atau "Langkah selanjutnya"

6. Memeriksa Tujuan Saat Ini

- Perintah : destinasi atau lokasi

- Penggunaan : Untuk memeriksa tujuan saat ini.

- Contoh : "Apa tujuan saya?" atau "Lokasi saat ini"

7. Berbagi Lokasi

- Perintah : bagikan

- Penggunaan : Untuk berbagi lokasi Anda dengan orang lain.

- Contoh : "Bagikan lokasi saya"

8. Mengganti Bahasa

- Perintah : bahasa atau ganti

- Penggunaan : Untuk mengganti bahasa antara bahasa Inggris dan Indonesia.

- Contoh : "Ganti bahasa" atau "Ubah ke bahasa Inggris"




Using Commands

To use the commands effectively:

- Activate the audio command feature by double-tapping the screen.
- Clearly speak the command and destination.
- Wait for the app to process and respond to your command.

Enjoy navigating with ease using our application! If you encounter any issues or need further assistance, please refer to the help section within the app or contact our support team.


On-Boarding Screen

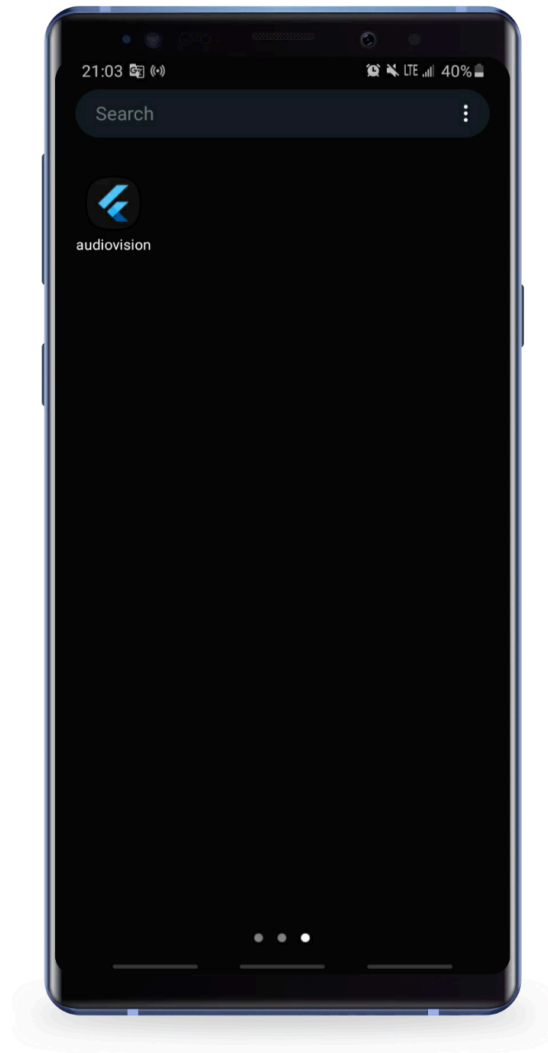

To get started, tap the "AudioVision" app on your smartphone.




Follow the on-screen guide to help you get set up quickly and easily.

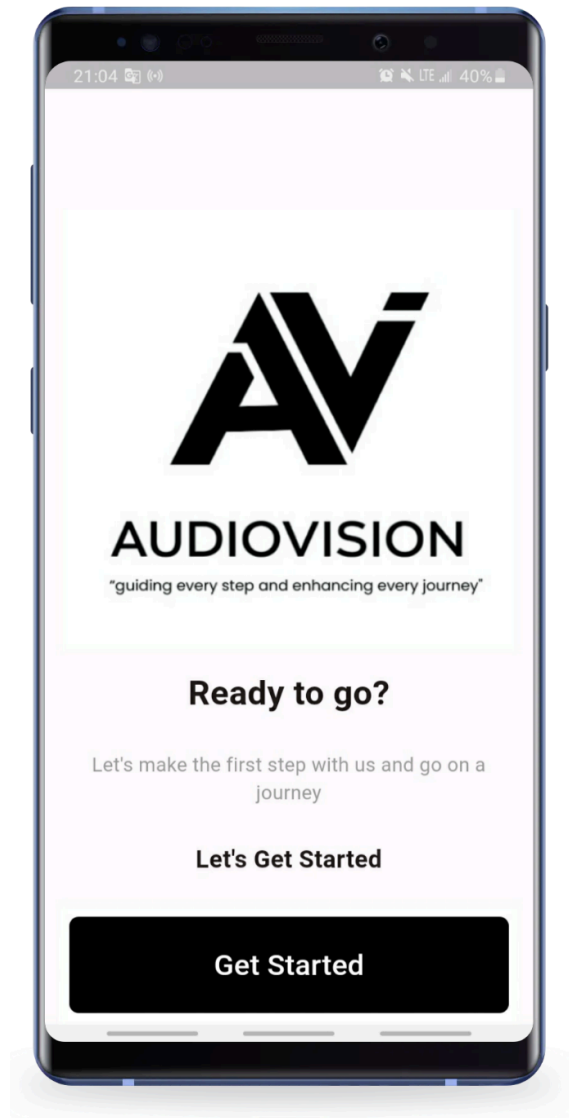


Once the Initial Guide Screen is completed, you'll hear instructions to click the "Get Started" button.

Home screen

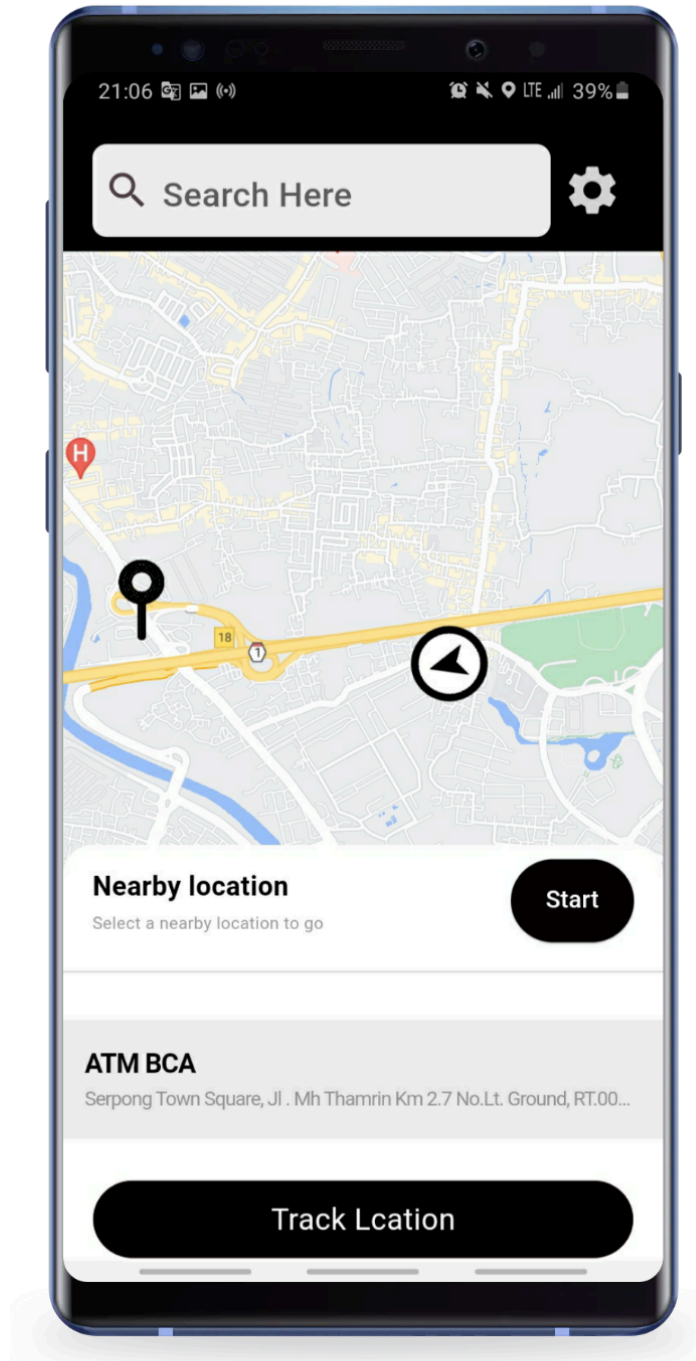


After clicking the "Get Started" button, you'll be directed to the Maps screen. Here, you'll see a map displaying your current location.


Start Navigate

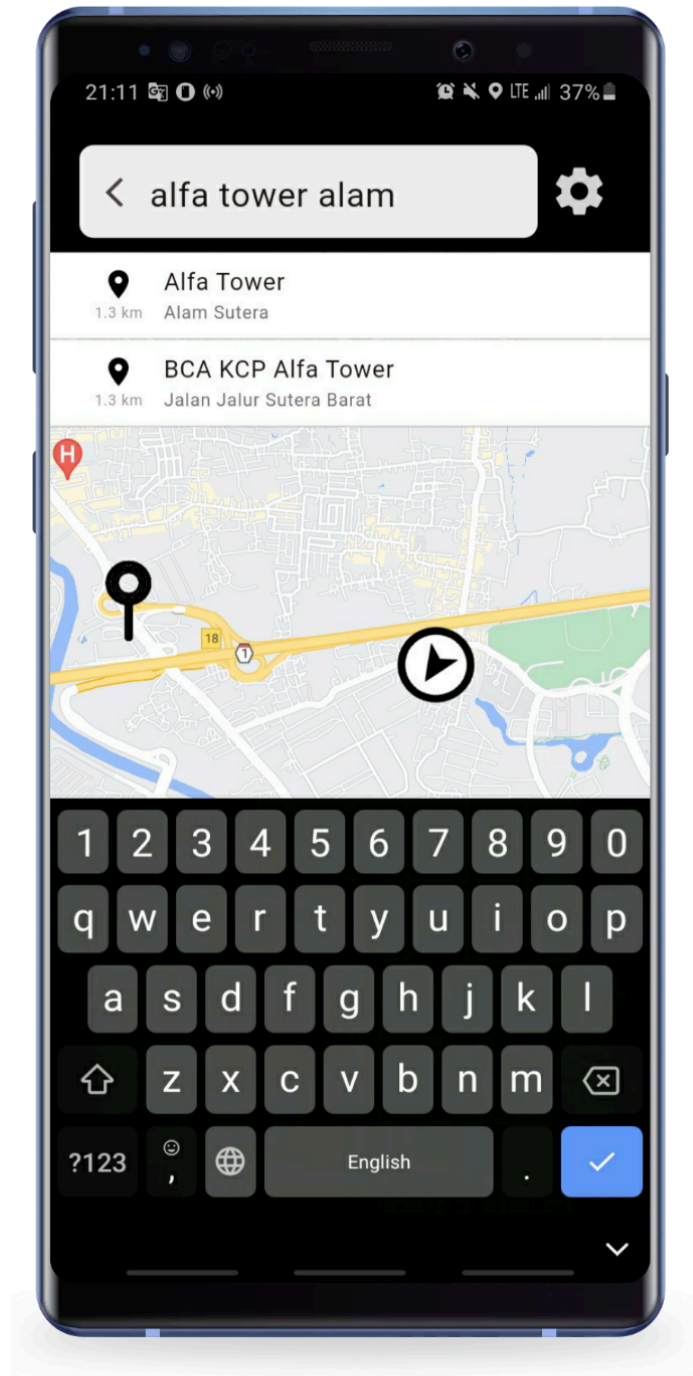

To start navigation, you need to input your destination first.

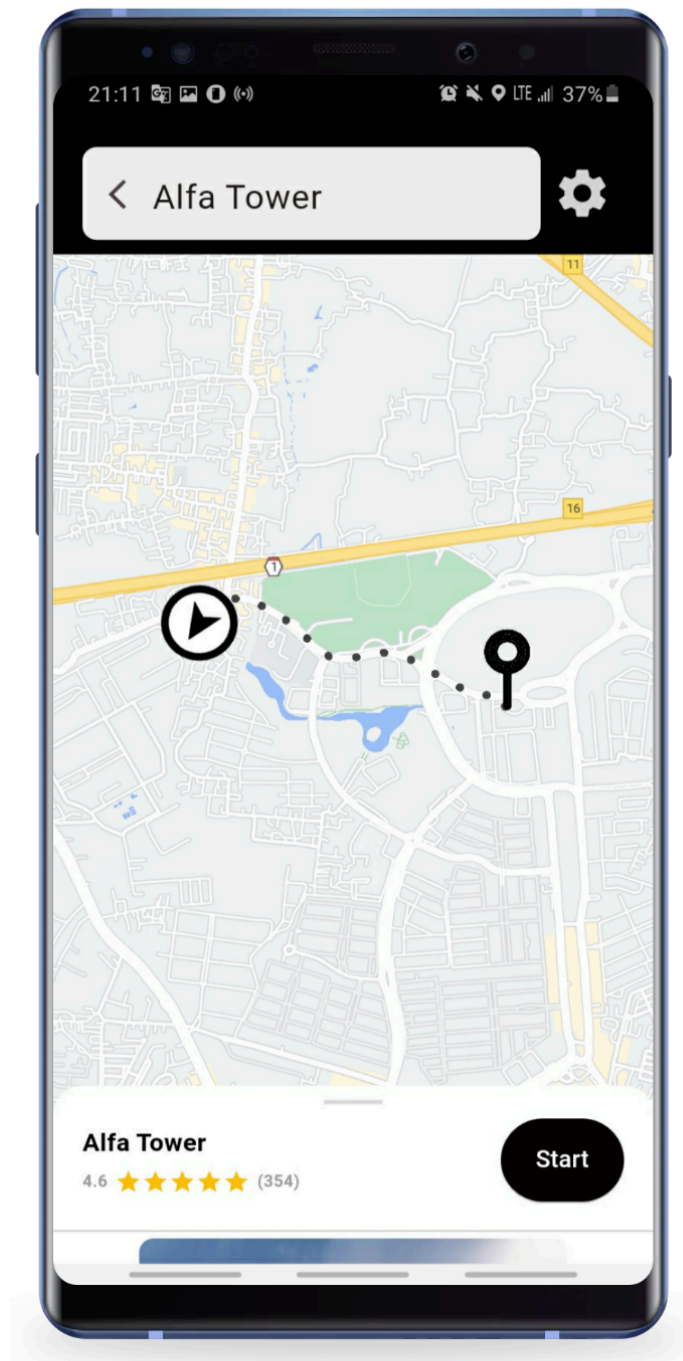


Once you've successfully entered your destination, the app will display a map with the route from your current location to your destination.
Don't worry, the route provided is tailored for pedestrians.
Now, you can start your journey by pressing the "Start" button.

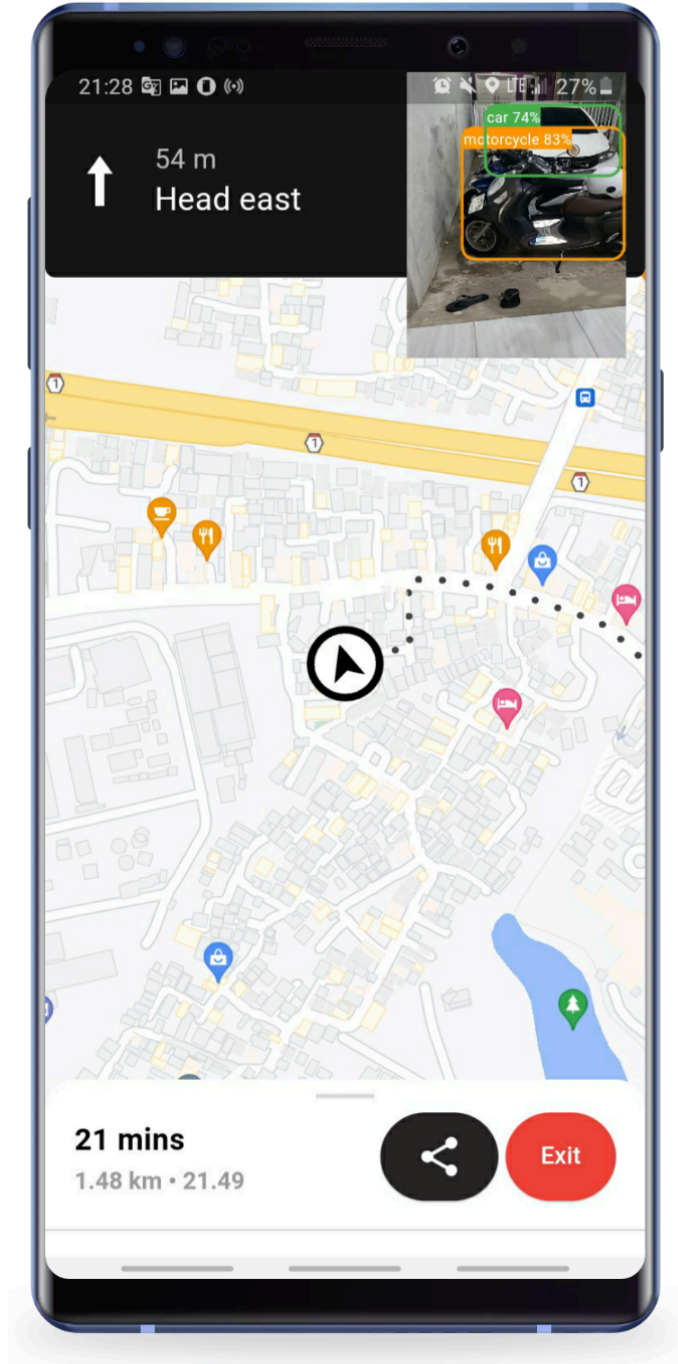


After you press the "Start" button, the AudioVision app will provide instructions about the route you'll take, and the camera will begin scanning for Real-time Object Detection.

At this stage, the app will provide information about detected objects. If a detected object is deemed hazardous, the app will issue a sound alert and vibration.


Share location

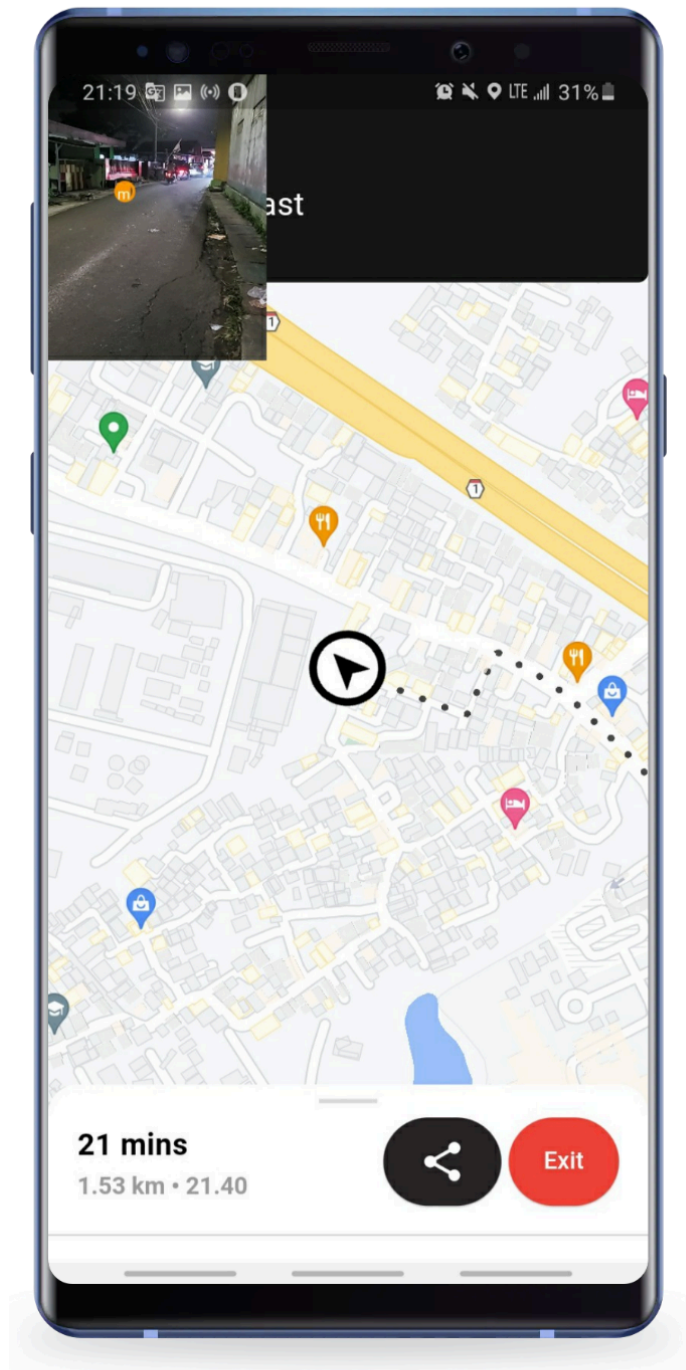


To use this feature, you simply need to press the "share" icon button.

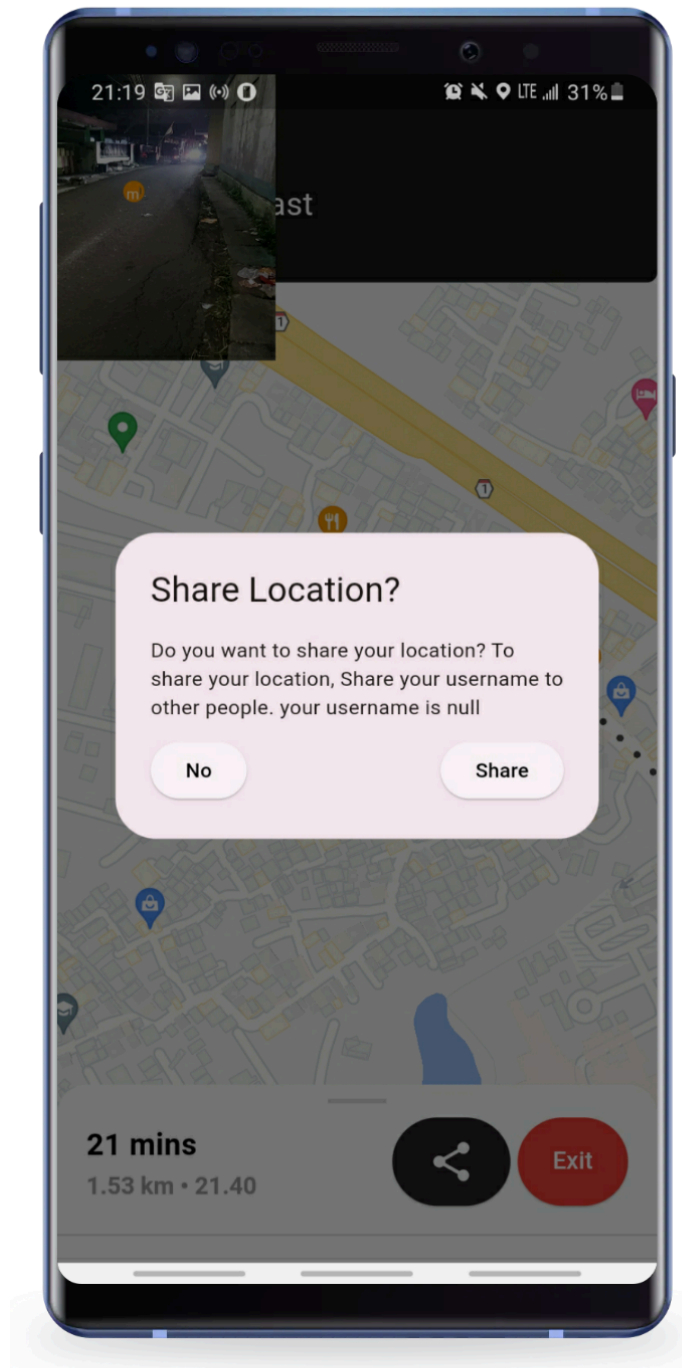

Afterward, a pop-up message will appear. Click the "share" button to share your location.

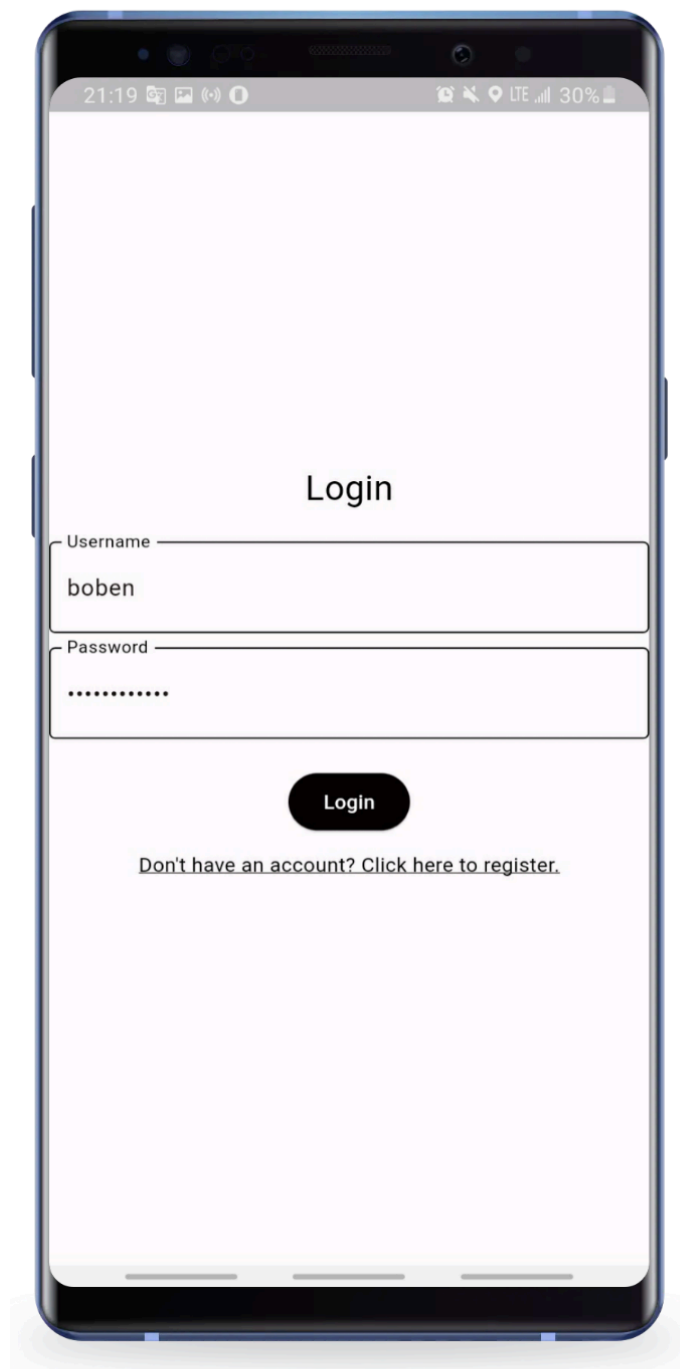


After that, you will be prompted to log in first. If you don't have an account yet, simply register first.
Click the "login" button.

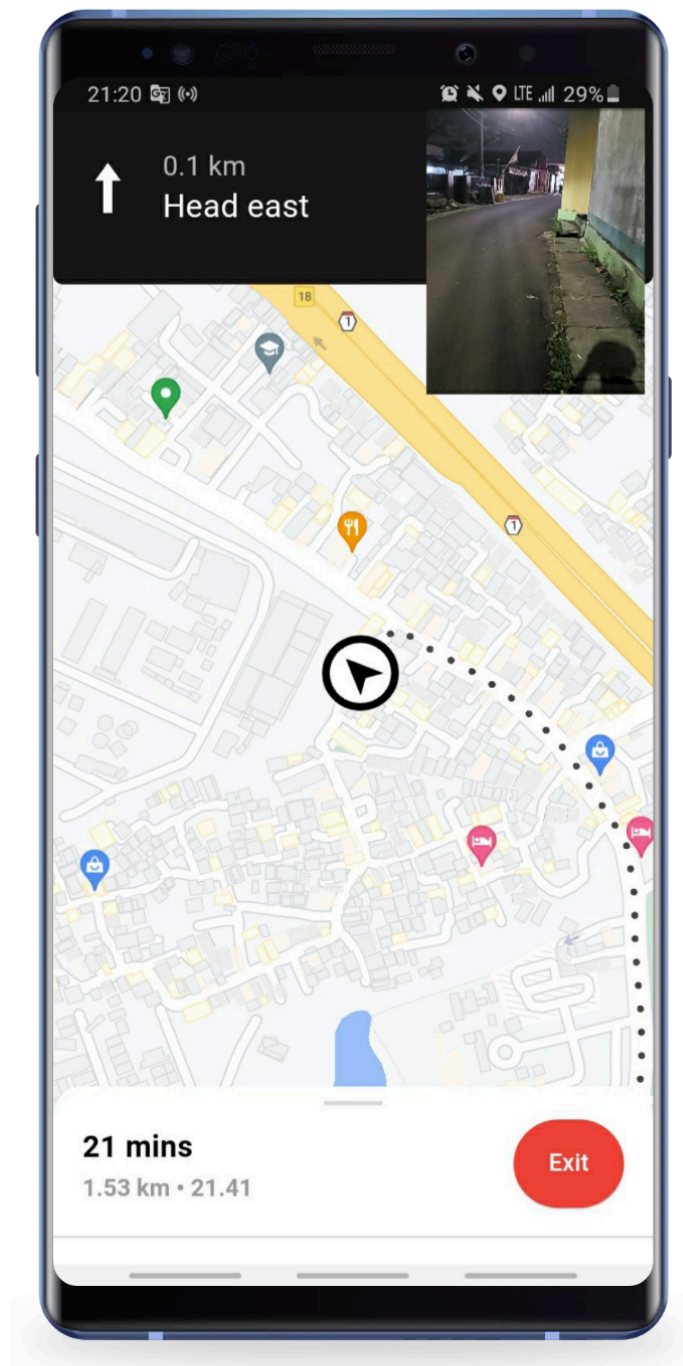


At this stage, you just need to share the username you registered during the registration process with trusted individuals as Observers. Once you reach your destination, your sharing location will stop.


Tracking Location

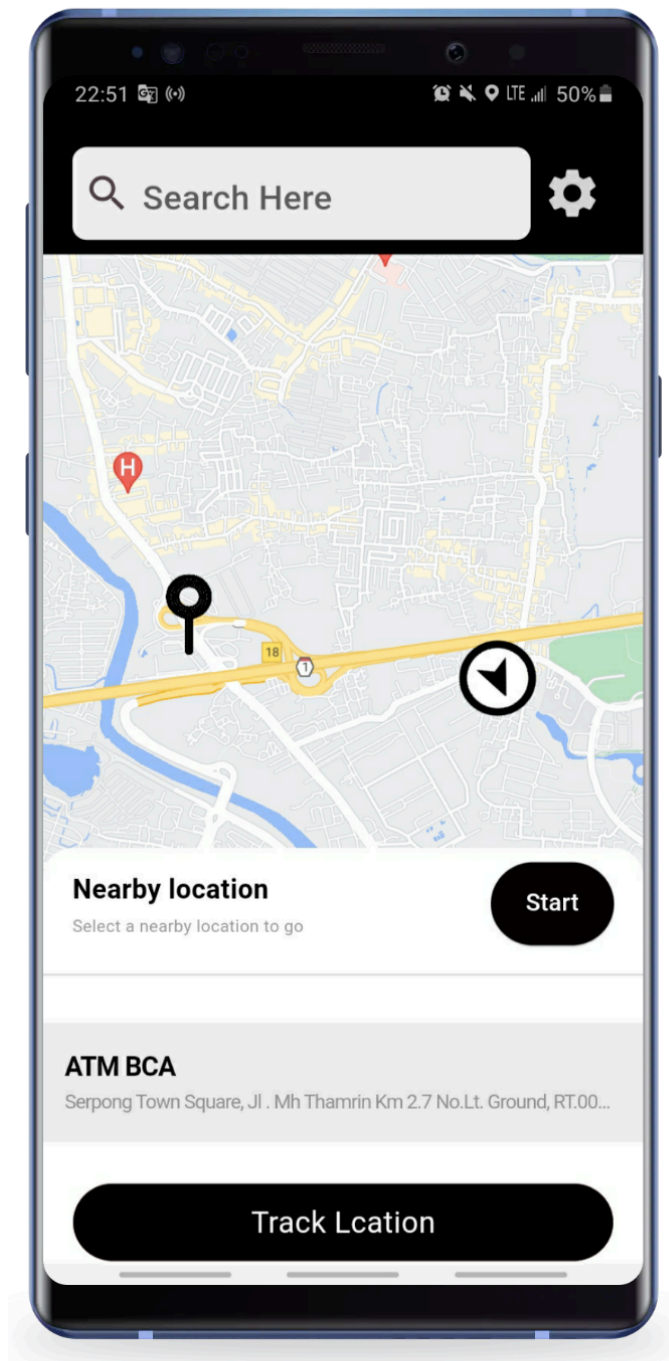


If you want to keep an eye on your homies, you can monitor their locations in real-time. Don't worry, you can do it easily. On the Homescreen, simply click the "Track Location" button.

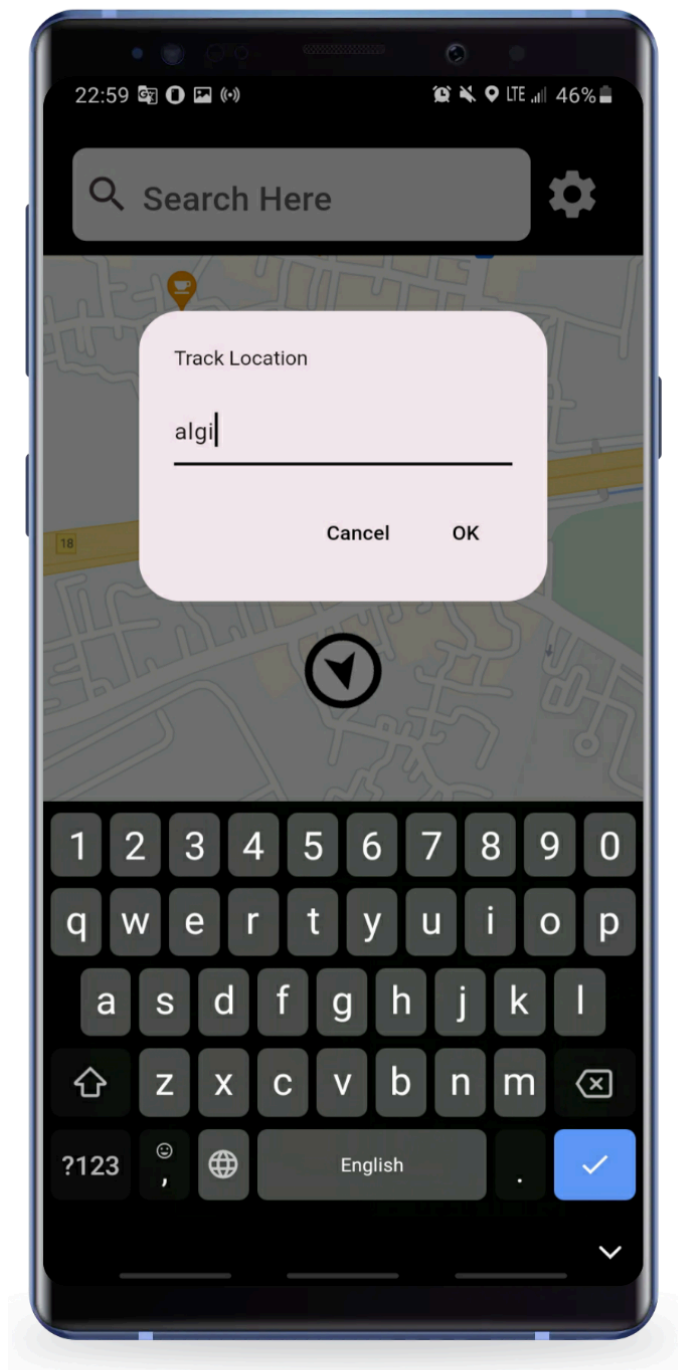

Next, you'll need to input the username used by your friend.

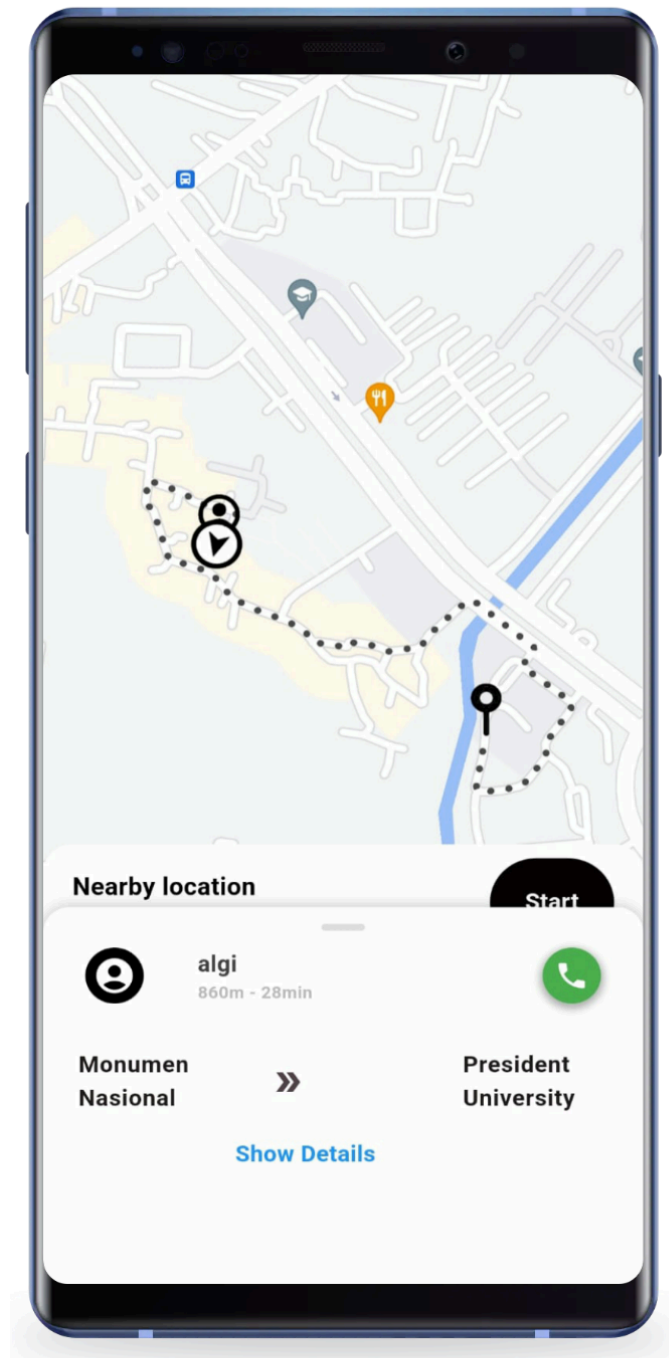


Now, you can track where your homies are heading and where they've been. You can even make a direct call to Emergency by tapping the "Call" icon button if something strange happens while tracking is taking place.

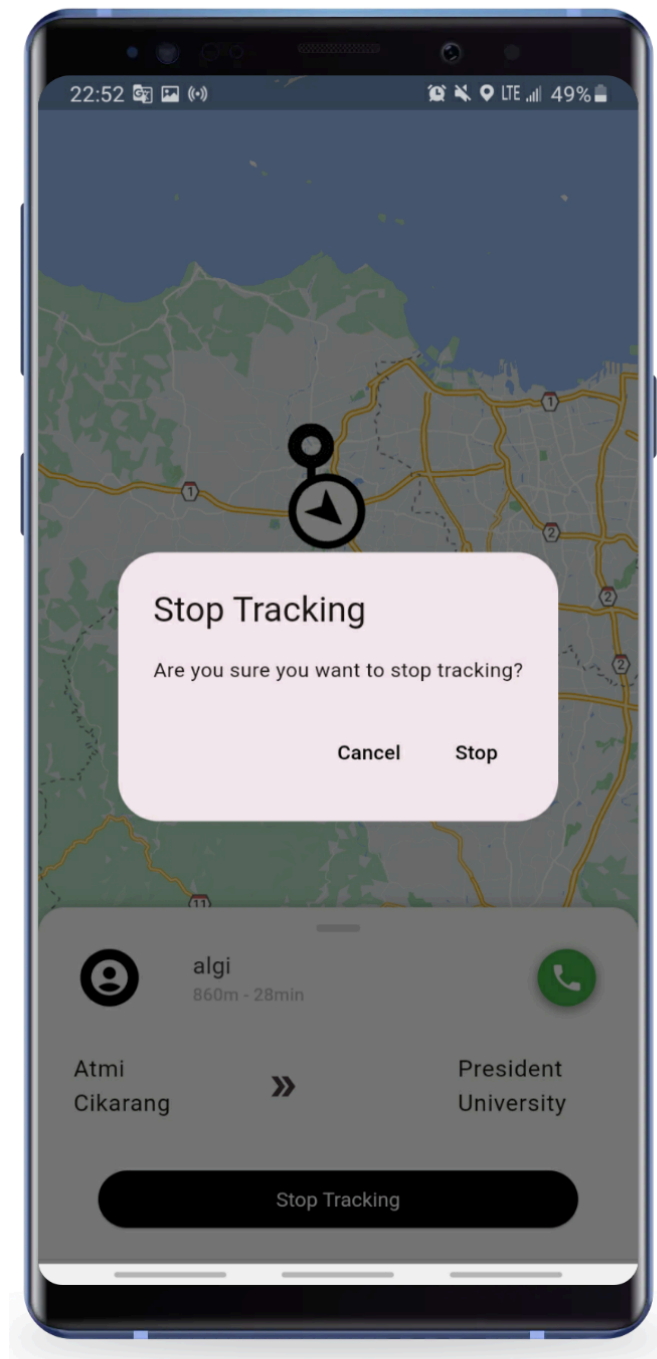


If you want to stop tracking their location, simply click the "Stop Tracking" button. Once the pop-up message appears, select "Stop". Your tracking session will successfully end.

Setting

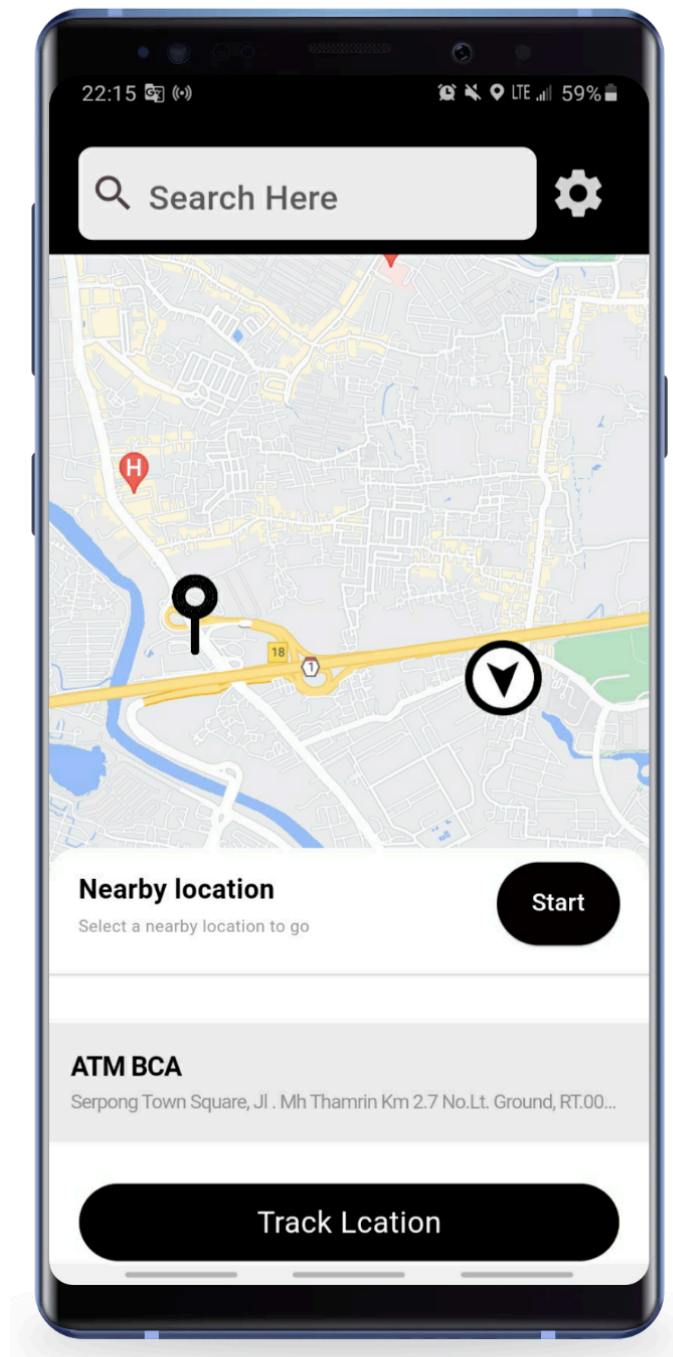

To access the Settings page, start by tapping on the "gear-shaped icon" located on the Homescreen

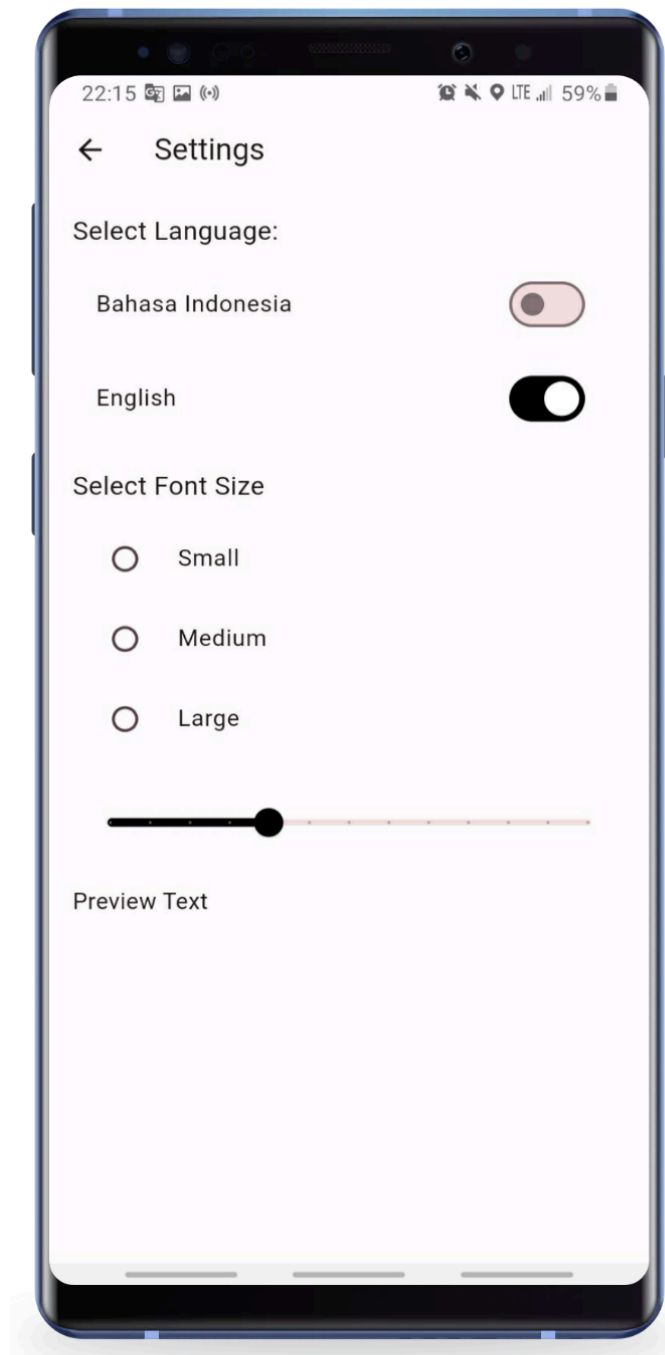

In the Settings page, customize your AudioVision experience by choosing your preferred language and adjusting the font size to match your preferences. It's all about making the app work for you!



Frequently Asked Questions (FAQ)

General Questions:

Q1: What is AudioVision?

A1: AudioVision is a revolutionary app designed to help visually impaired individuals navigate their daily activities safely and easily using advanced technology.

Q2: On which platforms is AudioVision available?

A2: For now!!! AudioVision is available only on Android platforms.

Installation and Setup:

Q3: What are the system requirements for installing AudioVision?

A3: You need a device with Android 8.0 (Oreo) or later. The device should have at least 4GB of RAM, 500MB of free storage space, a camera, a microphone, and an active internet connection.

Q4: How do I install AudioVision on my device?

A4: You can install AudioVision following steps on the “Installation Instructions” section. Follow the on-screen instructions to complete the installation and set up the Application.

Q5: I’m having trouble installing the app. What should I do?

A5: Ensure your device has enough storage space, a stable internet connection, and meets the minimum system requirements. Check app permissions in your device settings. If problems persist, don't hesitate to contact our support team.



Using AudioVision:

Q6: How do I use the Real-time Object Detection feature?

A6: Real-time object detection will be available when users have to “Start Navigate”. Open the app and Search Destination, then Start Navigate. Point your device's camera at your surroundings, and the app will provide audio feedback about the objects detected.

Q7: How do I give voice commands to the app?

A7: Simply activate the Audio Command feature and speak your command clearly. The app will respond and execute the command as instructed.

Q8: Can I track someone's location using AudioVision?

A8: Yes, with the Live Tracking feature, you can track the location of a user in real-time, provided they have shared their location with you.

Troubleshooting:

Q9: The app is not detecting objects correctly. What should I do?

A9: Ensure your camera lens is clean and unobstructed. Make sure there is enough lighting in your surroundings. If the issue persists, restart the app or your device.

Q10: I'm not receiving audio feedback from the app. How can I fix this?

A10: Check your device's volume settings and ensure the app has permission to use the microphone. Restart the app and try again. If the problem continues, contact our support team.



Contacting Support

If you encounter any issues or have questions while using AudioVision, our support team is here to help. Here's how you can reach us and find additional resources:

Contacting Support Team

Email Support

For any technical issues, account problems, or general inquiries, you can email our support team at:

- **Email:** AudioVision.helpdesk@outlook.com

Please include the following information in your email to help us assist you more efficiently:

- Your name
- Description of the issue
- Device type and operating system (e.g., Android 10,Ram 4 GB).
- Screenshots or error messages, if applicable

Feedback and Suggestions

We value your feedback and are constantly striving to improve AudioVision. If you have any suggestions for new features or improvements, please let us know.

- Don't hesitate to give feedback or suggestions, you can reach us via Our Emails.

Thank you for using AudioVision. We are committed to providing you with the best possible experience and support.

